


RPI Part #CME004 Replacement Filter Element INSTALLATION INSTRUCTIONS

- 1) Do not remove the Main Filter Housing when replacing the new Filter Element (RPI Part #CME004)
- 2) Press down on the two tabs to release the cap.
- 3) Inspect the seals located at each end. Replace the O-ring located in the End Cap.
- 4) Remove any debris before installing the new Element.
- 5) Peel off the sticker below. Place it on the Element to cover the bottom portion of the old sticker as shown in the illustration.

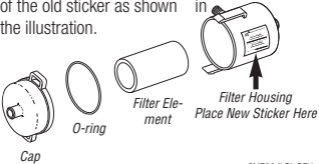

----- PEEL OFF STICKER -----

CME004LBL REV-

For Replacement Filter Element, order:
RPI Part #CME004
Replacement Filter Element

Date Installed: _____