



"The Alternate Source™"

replacement parts industries, inc.


CALL 800-221-9723 or 818-882-8611 • FAX 818-882-7028
E-MAIL: techsupport@rpiparts.com • WEBSITE: www.rpiparts.com

MIK078 STEAM TRAP KIT (New Type) INSTALLATION INSTRUCTIONS

Operation: This valve will seal off at the boiling point (212°F) and intermittently open again for approximately 5 seconds.

The Grommet (RPI Part #RXG009) is to be installed into the reservoir before the upper tube as shown below.

When servicing the **M9 UltraClave®**,
use enclosed **shorter** Lower Tube.

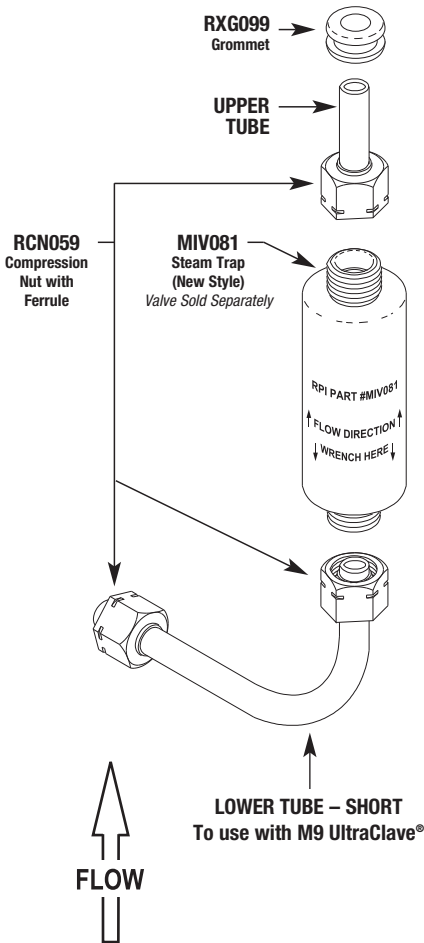

When servicing the **M11 UltraClave®**,
use enclosed **longer** Lower Tube.

