





Fits
**AUTOCRIT &
MICRO-HEMATOCRIT**

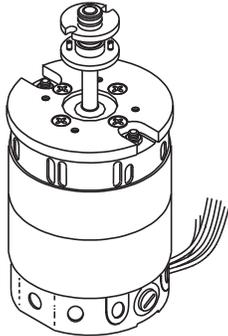

RPI Part #CAM115
OEM Part #42055801
(Old OEM Part #0551-600-000)

MOTOR & SPINDLE ASSEMBLY
• 120VAC, 60 Hz.
• 11,700 RPM (+22%, -10%)


FREE!
Just ask for it!

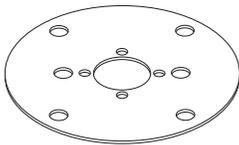

MOTOR MOUNTING PLATE
RPI Part #CAM049B

Fits: Motor & Spindle Assembly


Fits
**AUTOCRIT II &
MICRO-HEMATOCRIT II**

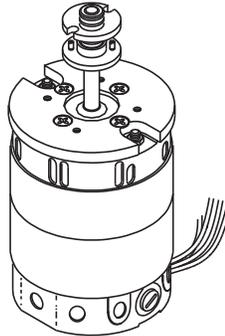

RPI Part #CAM115
OEM Part #42055801
(Old OEM Part #0558-600-00)

MOTOR & SPINDLE ASSEMBLY
• 120VAC, 60 Hz.
• 11,700 RPM (+22%, -10%)


FREE!
Just ask for it!

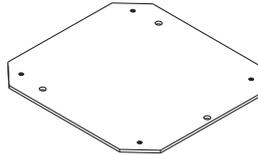

LOWER CHASSIS COVER PLATE
RPI Part #CAM115B

Fits: Chassis


Fits
**AUTOCRIT ULTRA III &
MICRO-HEMATOCRIT ULTRA III**

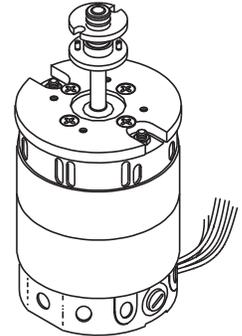

RPI Part #CAM115
OEM Part #42055801
(Old OEM Part #0575-600-000)

MOTOR & SPINDLE ASSEMBLY
• 120VAC, 60 Hz.
• 11,700 RPM (+22%, -10%)


FREE!
Just ask for it!

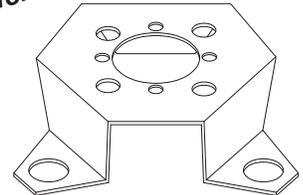

MOTOR MOUNTING PLATE
RPI Part #CAM050B

Fits: Motor & Spindle Assembly

FREE!
Just ask for it!

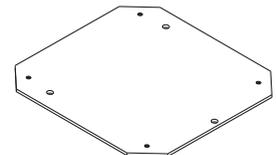

LOWER CHASSIS COVER PLATE
RPI Part #CAM115B

Fits: Chassis

Important Note Regarding Autocrit and Micro-Hematocrit: The head adapter on the old style RPI motors that fit the Autocrit and Micro-Hematocrit can not be removed to allow reuse of the mounting plate. So when replacing the old style RPI Motor Assembly (RPI Part #CAM049) with the new style (RPI Part #CAM115), you will need to replace the round style mounting plate. RPI provides a new round style mounting plate free of charge when you order the the new style Motor Assembly (RPI Part #CAM115).

Important Note Regarding the Autocrit II and Micro-Hematocrit II: The head adapter on the old style RPI motors that fit the Autocrit II and Micro-Hematocrit II can not be removed to allow reuse of the mounting plate. So when replacing the old style RPI Motor Assembly (RPI Part #CAM050) with the new style (RPI Part #CAM115), you will need to replace the lower chassis cover plate. RPI provides a new lower chassis cover plate free of charge when you order the the new style Motor Assembly (RPI Part #CAM115).

Important Note Regarding the Autocrit III and Micro-Hematocrit III: You will need to replace the triangular style mounting plate and a lower chassis cover plate. RPI provides new plates free of charge when you order the new style Motor Assembly (RPI Part #CAM115).